



Confidence in posture solutions

# Azalea<sup>®</sup> Max



A tailor-made solution for heavier individuals



Yes, you can.<sup>®</sup>

The **Rea Azalea Max** is designed to meet the specific needs of people with extra weight. For effortless positioning we recommend power tilt and backrest recline. The **Rea Azalea Max** is a niche product which complements the existing **Azalea** family from Invacare.

## Continued Compatibility

The **Rea Azalea Max** includes features such as tilt-in-space with weight-shifting and DSS (Dual Stability System) mechanism from the **Azalea** range. Stability in this wheelchair originates from body length compensation in the back and legrests, beneficial for those with limited movement who have difficulty changing position.

## Features and Options



### Two-way backrest

The lower and upper parts of the shaped backrest are reversible, giving the option of more space in the seat or chest area.



### Footplate converter

Divided into three parts, the footplate comes with a removable mid section, to facilitate transfers.



### Power tilt and recline

With the hand control the person or the attendant can easily change seating position.



### High ergonomic push handle

The ergonomically shaped push handles are positioned to optimise the pushing for the attendant. An alternate version with higher handles offers a different position for the hands and is available as an option.

## Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk)

	Seat width	Seat depth	Seat height	Backrest height	Armrest height	Length of legrests	Seat angle	Backrest angle
<b>Rea Azalea Max</b>	550; 610 + 50 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 500 mm	1° - 20°	0°-30°
<b>Rea Azalea Max Assist</b>	550; 610 + 100 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 500 mm	1° - 20°	0°-30°

	Driving unit width	Total height	Total length inc. Footrests	Total weight	Max. User weight	Transport weight	Maximum slope for use of parking brake
<b>Rea Azalea Max</b>	SW + 250 mm	1020 - 1240 mm	1050 - 1120 mm	54 kg	max. 180 kg	33 kg	7°
<b>Rea Azalea Max Assist</b>	SW + 220 mm	1020 - 1240 mm	1050 - 1120 mm	54 kg 1	max. 180 kg	37 kg	7°

1. Transport weight increases by 4,2 kg if Power tilt & recline is selected.

## Frame colour

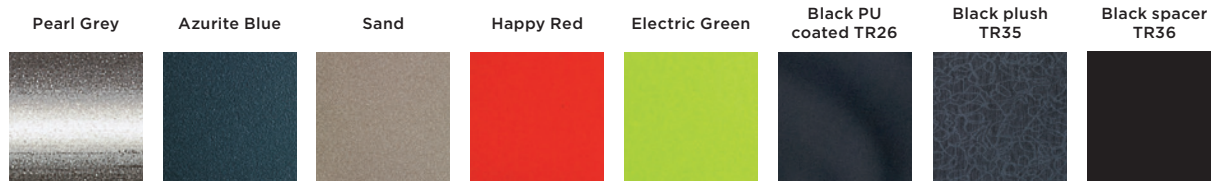

## Upholstery



Invacare Limited  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail: [ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

