

rea[®]

Clematis[®] Pro



Supreme comfort | Effortless adjustability | Compact design



An evolution in tilt-in-space passive wheelchairs

Our new Rea Clematis Pro has been designed to enhance the experience for both client and care giver. It's perfectly suited for those individuals who need the support and comfort of a tilt in space wheelchair. Part of our market leading Rea Passive range, the Clematis Pro offers a high level of postural support, with 25° of seat tilt and 30° of backrest recline for effortless repositioning. A wide range of accessories helps deliver exceptional comfort, support and stability for those who spend large amounts of time in their wheelchair.

Supreme comfort backrest

The rigid plastic backrest shell is height adjustable and incorporates 30° of backrest recline offering a high level of support where it's needed. The addition of the Laguna backrest cushion with added lateral support, ensures a comfortable and stable sitting position.

DSS® (Dual Stability System)

The unique weight shifting technology DSS, ensures excellent stability of the wheelchair. When tilted and reclined, the seat moves forward keeping the clients weight in the footprint of the chassis, ensuring maximum safety and increased confidence.

Compact design

The compact frame design offers superb manoeuvrability, especially indoors. The rear wheels have been moved forward to allow better access to the handrims for those who can propel themselves.



7 FEATURES



1 Adjustable backrest

- ▶ Quick and easy to set up and adjust to individual needs
- ▶ Backrest cushion with lateral support for added comfort
- ▶ Precise fit for a range of postural requirements

2 Larger seat tilt

- ▶ Max seat tilt of 25°
- ▶ Ensures improved pressure redistribution
- ▶ Reduces shear forces

3 Adjustable seat plate

- ▶ Adjustable in depth

4 Angle adjustable legrests

- ▶ Allows the legrest to be elevated for comfort when in a tilt and/or recline position
- ▶ Height and depth adjustable calfpad
- ▶ Height, depth, angle and width adjustable footplate to find the perfect position for the feet

5 Compact frame

- ▶ Easy to manoeuvre in the tightest of spaces
- ▶ Optimised for indoor use
- ▶ Anti-tippers to improve safety

6 Dual Stability System

- ▶ Enhanced stability in every tilt and recline position
- ▶ High pivot point to minimise shear forces.

7 Width adjustment

- ▶ The double armrest receiver offers width adjustment of 5 cm as standard
- ▶ Increased space in the seat area

Additional features

- ▶ **New concave armpad**
 - The one piece moulded concave armpad ensures the arm stays in position throughout the tilt and recline movements.
 - Easy to clean for better infection control
- ▶ **New style Laguna backrest cushion**
 - Improved comfort and postural support



Range of other options include

- Table tray
- Transportation kit
- Removable pommel
- Adjustable Trunk supports
- Tilt lock mechanism
- Foot operated brake (Transit only)
- Headrest cover
- Footplate cover
- Spokeguards

Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.



Seat width ▼

390 - 540 mm
in increments of
50mm




Seat depth ▼

430 - 500 mm



Seat height ▼

400 - 450 mm

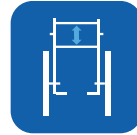

Backrest height ▼

570 mm -710 mm



Armrest height ▼

230 - 350 mm

Rea Clematis Pro



Length of footrests ▼

355 - 465 mm




Seat angle ▼

0° - 25°



Backrest angle ▼

0° - 30°

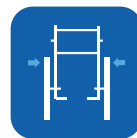

Overall width ▼

SW + 190 mm



Total height ▼

1200 - 1550mm

Rea Clematis Pro



Total length inc.
footrests ▼

995 mm - 1470 mm



Total weight ▼

from 30 kg



Max. user weight ▼

135kg



Transport weight ▼

from 20 kg




Maximum slope
for use of parking
brake ▼

7°

Rea Clematis Pro

Frame colour ▼

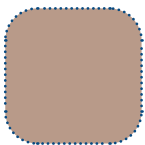

**Pearl
Beige**



Black TR26

Upholstery ▼



Dual Stability System

Rea, Rea design and DSS (Dual Stability System) design are registered trademarks of Invacare International. The Rea Clematis Pro is designed with DSS (Dual Stability system). It provides the user with enhanced stability and safety when being tilted or the backrest is reclined.

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie



Yes, you can.®

© 2019 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Rea Clematis Pro - EU - 08/2019 - Part number 1632922

